+43°C

12 RESPONSIBLE
CONSUMPTION
AND PRODUCTION

Vestfrost Solutions is working towards reaching the UN - Global Sustainable Development Goals by 2030.

The Sustainable Development Goals are the blueprint to achieve a better and more sustainable future for all. In order to implement Goal no 12 "Responsible Consumption and Production", this manual has been printed on recycled paper.

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WARNING

As the appliance contains flammable refrigerant, as stated on nameplate, it is essential to ensure that the refrigerant pipes are not damaged.

The quantity and type of the refrigerant used in your appliance is indicated on the rating plate.

Standard EN378 specifies that the room in which you install your appliance must have a volume of 1m³ per 8 g of hydrocarbon refrigerant used in the appliances. This is to avoid the formation of flammable gas/air mixtures in the room where the appliance is located in the event of a leak in the refrigerant circuit.

WARNING:

Keep ventilation openings in the appliance or in built-in structures must be kept clear.

WARNING:

Do not use other mechanical devices or means to accelerate the defrosting process or to remove rime other than those recommended by the manufacturer.

WARNING:

Do not damage the refrigerant system.

WARNING:

The appliance may not be exposed to rain.

WARNING:

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack experience and knowledge, unless they have been given supervision or instructions concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.

WARNING:

Children must not play with, on, or around the appliance.

WARNING:

Children must not clean the appliance or carry out general maintenance unless they are at least 8 years old and are being supervised.

WARNING:

Always, keep the keys in a separate place and out of reach of children.

WARNING:

Before servicing or cleaning the appliance, switch off circuit breaker.

**WARNING:**

Danger risk of fire or explosion. Flammable refrigerant used, as stated on nameplate. To be repaired only by trained personnel.

WARNING:

Do not store explosive substances such as aerosol cans with a flammable propellant in this appliance.

WARNING:

When positioning the appliance, ensure the power cord is not trapped or damaged.

WARNING:

Do not locate multiple portable socket-outlets or portable power supplies at the rear of the appliance.

WARNING:

Appliance use flammable insulation blowing gas. For information about safe disposal, please contact your local disposal service. See section for Disposal.

Contents

WARNING	2
Periodic preventive maintenance checks ...	5
Spare part list MF 114	6
Spare part list MF 214	8
Spare part list MF 314	10
Vital components	12
Health and safety guidance – Warning!....	13
Required basic tools	14
Controller	15
Temperature adjustment.....	15
Starting device replacement.....	17
Thermostat sensor replacement.....	20
Thermostat replacement.....	22
Thermometer replacement	24
Adjustment of lid	27
Wire diagram	28
Compressor replacement	29
On-site checklist	30
Trouble shooting	31
Technical support.....	32
Recycling procedures	33

Periodic preventive maintenance checks

Daily Check:

Monitor Temperature.
Lid fits and lock tight to cabinet
Lid gasket not faulty.

Weekly maintenance:

Check for the ice layer
Defrost equipment if ice > 1 cm
Clean the lid.

Monthly maintenance:

Clean grille for compressor compartment.
Clean the refrigerator with lukewarm water
and mild detergent.


Yearly maintenance:

Check electrical connections and components.

PQS Code	Model	PQS Performance specifications Specification reference:	PQS Independent type-testing protocol Product verification protocol:
E003/024	MF 114	E003/FZ01.2	E003/FZ01-VP.2
E003/025	MF 214	E003/FZ01.2	E003/FZ01-VP.2
E003/023	MF 314	E003/FZ01.2	E003/FZ01-VP.2


Spare part list MF 114

Position	Item number	Item name
0061	0-6538001	Filter drier,
0070	A921001	Cover + wirring clamp
0071	A921115	Starting device, 103N0021
0087	0-6038175	Base plate fittings, complete
0129	8-036038255	Compressor - NLE 9KTK , Complete, R600a
0550	6520317	Run cap. 117-7119 4uF/4,8 receptacles
0899	7020983	Sensor
3125	3010457	Pipe holder Ø6
5000	3010049	Drain plug
5002	6000959	Foot with pin
5057	6010084	WOT condenser
5108	5000917670	Lid foamed without handle/hinges
5375	6000965	Partition, complete
5376	5020105	Lid for fast freeze compartment
5377	3020245-23	Handle
5401	1510133	Hinge
5402	1510136	Hinge with spring
5450	3011135-01	Top part for hinge cover
5451	3010032-01	Bottom part for hinge cover
5605	304090501	Handle with lock (push and turn)
5606	8090342-94	Inlay for handle,
5682	1510046	Key, set - 2pc.
5683	2040145	Catch for handle
5684	3010265-01	Cover for catch
5685	600098801	Lock with keys,(push and turn)
5707	7020248	Voltage stabilizer
5710	7020239	Thermostat, freezer
5713	7020960	Sensor NTC S1, Cool
5738	7010027	Pilot lamp, green
5756	7020164	Thermometer, solar
5771	3010051-01	Thermostat knob
5783	3040146	Control panel
5870	6000308	Fan
5901	353090001	Basket with handle, white
5902	3010107-01	Spare handle for basket
5970	3010308-01	Motor screen
6553	7535360	Mains lead. Schuko. 3G1/2600

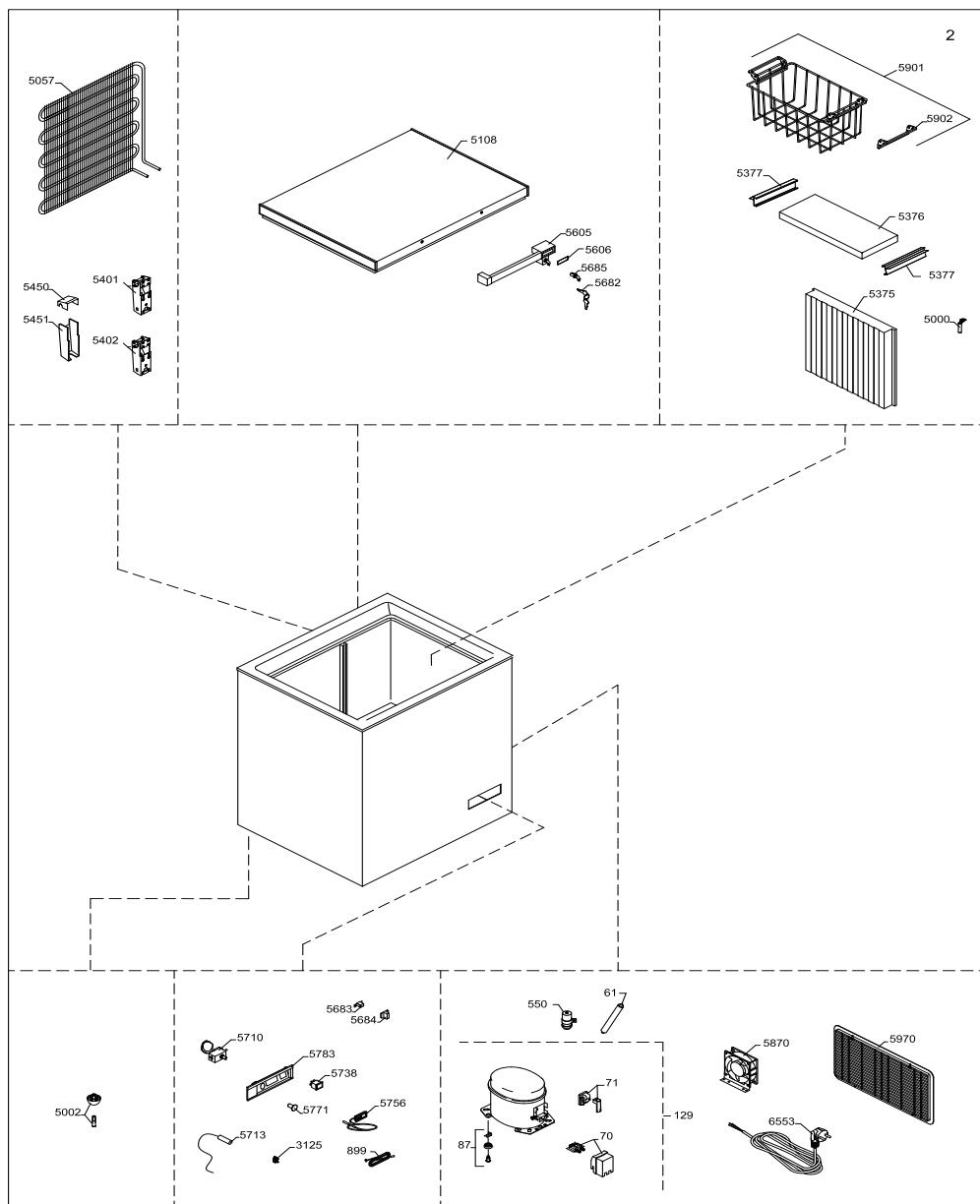
Spare part list MF 214

Position	Item number	Item name
0061	0-6538001	Filter drier,
0070	A921001	Cover + wirring clamp
0071	A921115	Starting device, 103N0021
0087	0-6038175	Base plate fittings, complete
0129	8-036038255	Compressor - NLE 9KTK , Complete, R600a
0550	6520317	Run cap. 117-7119 4uF/4,8 receptacles
0899	7020983	Sensor
3125	3010457	Pipe holder Ø6
5000	3010049	Drain plug
5002	6000959	Foot with pin
5057	6010084	WOT condenser
5116	5000918651	Lid foamed without handle/hinges
5374	3040184	Ice pack tray
5375	6000965	Partition, complete
5376	5020105	Lid for fast freeze compartment
5377	3020245-23	Handle
5401	1510135	Hinge with spring
5450	3011135-01	Top part for hinge cover
5451	3010032-01	Bottom part for hinge cover
5605	304090501	Handle with lock (push and turn)
5606	8090342-94	Inlay for handle,
5682	1510046	Key, set - 2pc.
5683	2040145	Catch for handle
5684	3010265-01	Cover for catch
5685	600098801	Lock with keys,(push and turn)
5707	7020248	Voltage stabilizer
5710	7020239	Thermostat, freezer
5713	7020960	Sensor NTC S1, Cool
5738	7010027	Pilot lamp, green
5756	7020164	Thermometer, solar
5771	3010051-01	Thermostat knob
5783	3040146	Control panel
5870	6000308	Fan
5901	353090001	Basket with handle, white
5902	3010107-01	Spare handle for basket
5970	3010308-01	Motor screen
6553	7535360	Mains lead, Schuko, 3G1/2600

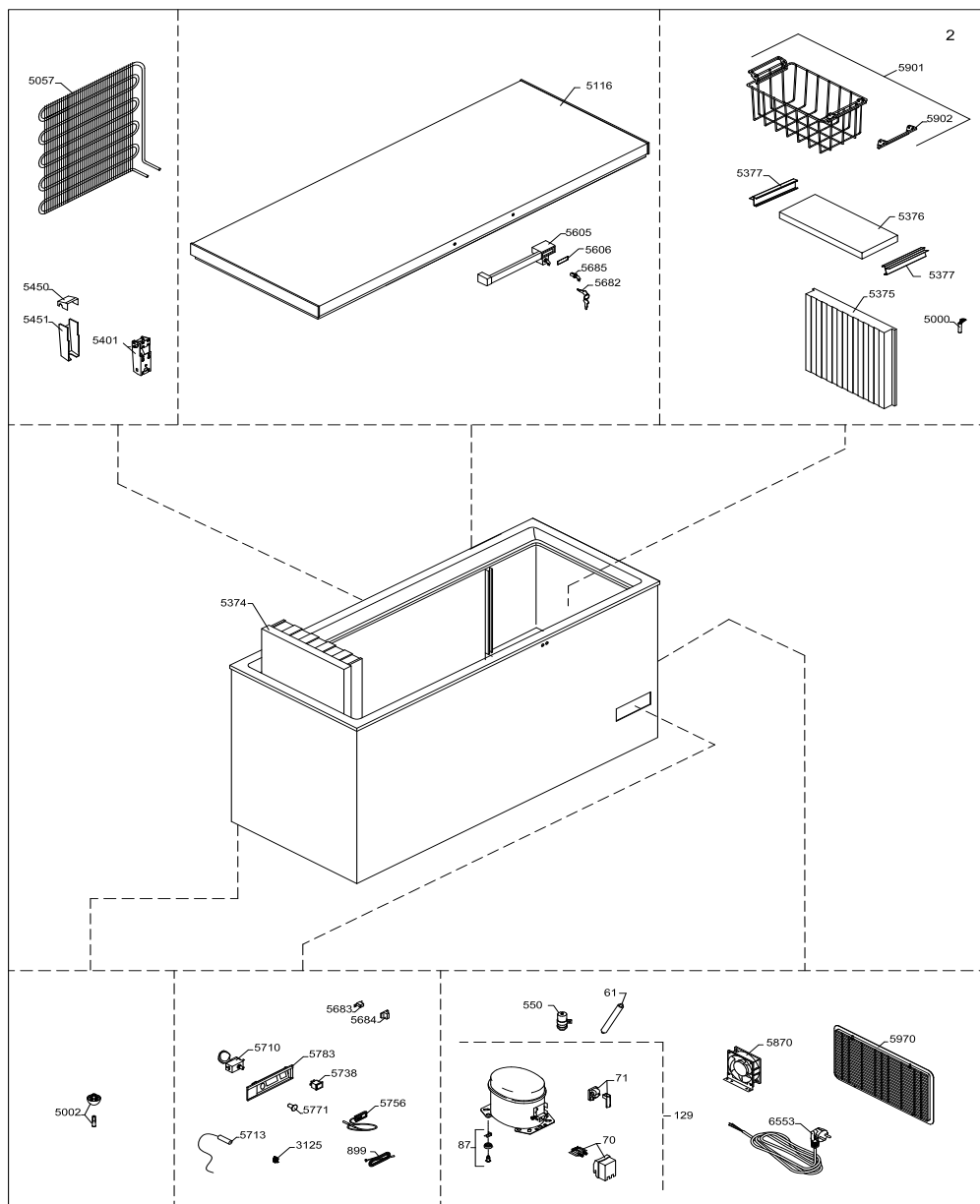Spare part list MF 314

Position	Item number	Item name
0061	0-6538001	Filter drier,
0087	0-6038175	Base plate fittings, complete
0129	8-03606510309	Compressor - HXK12AT , Complete,
0550	6520510	Run cap. 4uF/4,8 receptacles
0899	7020983	Sensor
1126	0-A921012	Cover + wirring clamp
1141	6520511	Starting device PTC
3125	3010457	Pipe holder Ø6
5000	3010049	Drain plug
5002	6000959	Foot with pin
5057	6010084	WOT condenser
5122	5000914651	Lid foamed without handle/hinges
5374	3040184	Ice pack tray
5375	6000965	Partition, complete
5376	5020105	Lid for fast freeze compartment
5377	3020245-23	Handle
5402	1510136	Hinge with spring
5450	3011135-01	Top part for hinge cover
5451	3010032-01	Bottom part for hinge cover
5605	304090501	Handle with lock (push and turn)
5606	8090342-94	Inlay for handle,
5682	1510046	Key, set - 2pc.
5683	2040145	Catch for handle
5684	3010265-01	Cover for catch
5685	600098801	Lock with keys,(push and turn)
5707	7020248	Voltage stabilizer
5710	7020239	Thermostat, freezer
5713	7020960	Sensor NTC S1, Cool
5738	7010027	Pilot lamp, green
5756	7020164	Thermometer, solar
5771	3010051-01	Thermostat knob
5783	3040146	Control panel
5870	6000308	Fan
5901	353090001	Basket with handle, white
5902	3010107-01	Spare handle for basket
5970	3010308-01	Motor screen
6553	7535385	Mains lead

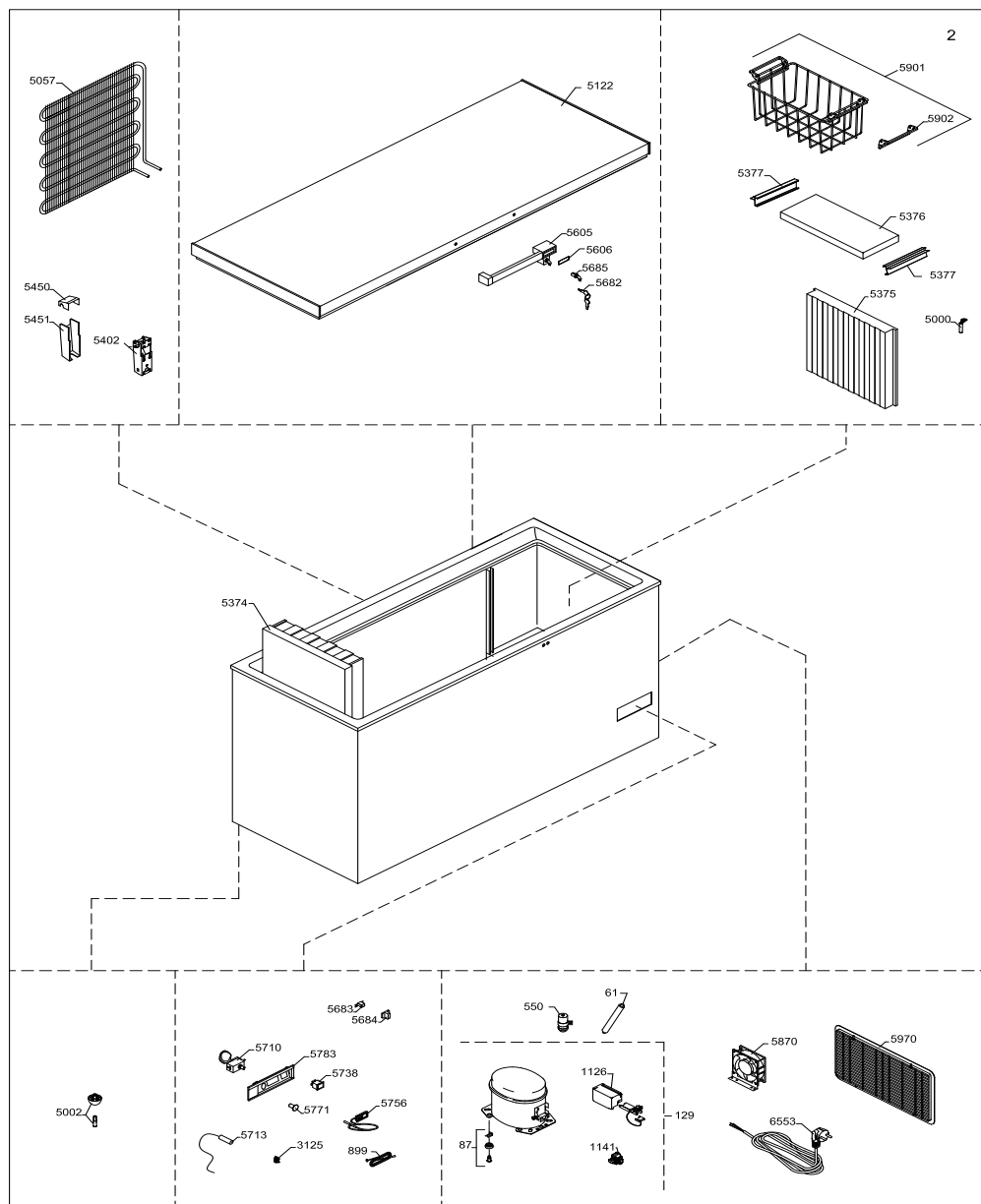Vital components


Position	Item no	Description
----------	---------	-------------

MF 114 & MF 214

0129	8-036038255	Compressor
------	-------------	------------


MF 314

0129	8-03606510309	Compressor
------	---------------	------------

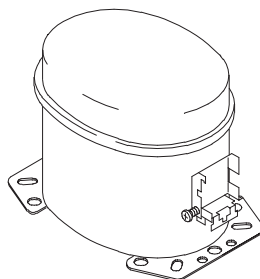
5713	7020960	Sensor
------	---------	--------

5710	7020239	Thermostat
------	---------	------------

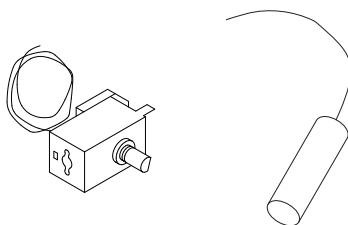
0550	6520317	Run capacitor
------	---------	---------------

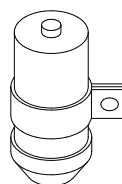MF 114 & MF 214

0070	A921001	Cover & wiring clamp
------	---------	----------------------


0071	A921115	Starting device
------	---------	-----------------

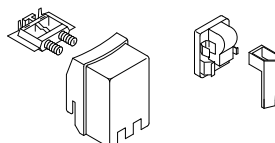MF 314

1126	0-A921012	Cover & wiring clamp
------	-----------	----------------------

1141	6520511	Starting device
------	---------	-----------------

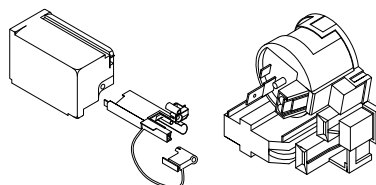Health and safety guidance – Warning!

Before any repair job be aware of following!

WARNING:

Before servicing or cleaning the appliance, disconnect it from power source.

WARNING:

Danger risk of fire or explosion. Flammable refrigerant used. To be repaired only by trained personnel.

Required basic tools

1. Flexible socket wrench

2. Socket wrench - size 7mm

3. Nose plier

4. Phillips screwdriver

5. Torx screwdriver - size t10+t20

6. Screwdriver - size 1,0x6,0

7. Screwdriver – size 0,6x3,5

8. Multimeter


Proposed additional service kit/items

Sealing kit

Tar tape

Extra self-tapping screws

Controller

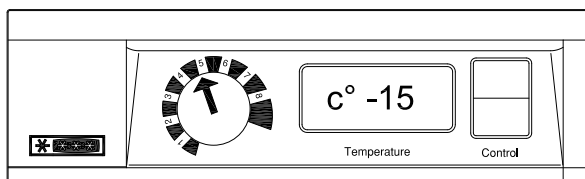
Before the appliance is loaded with vaccines ensure that the temperature is below -15°C :

- Set the thermostat at position 5
- Check the temperature on the thermometer
- If the temperature in the vaccine compartment is above -15°C the thermostat should be adjusted. The temperature in the vaccine compartment must always be monitored on the thermometer and be within the range -15 to -25°C .

Due to the mechanical tolerance of the thermostats you should always control the temperature during cool down .

Temperature adjustment



This appliance is equipped with one thermostat.

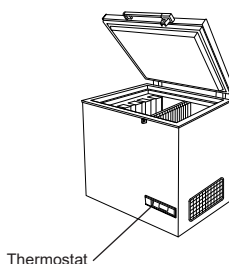
The thermostat is a electronic thermostat that controls the temperature in the vaccine and icepack freezing compartments. Setting 1 applies minimum cold and setting 8 maximum cold.

Motor compartment

How to get access to the motor compartment.

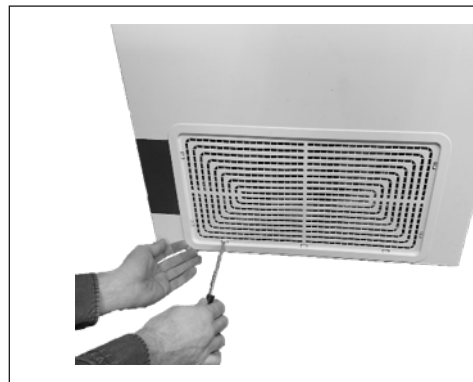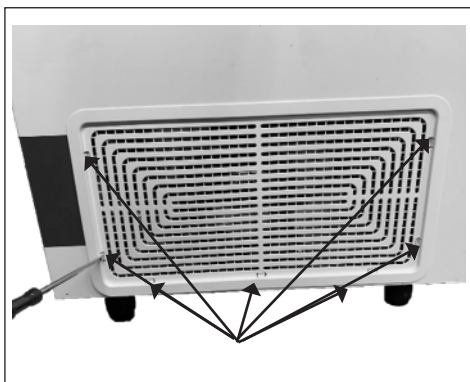1. Use a screwdriver to unlock all 7 clamps.

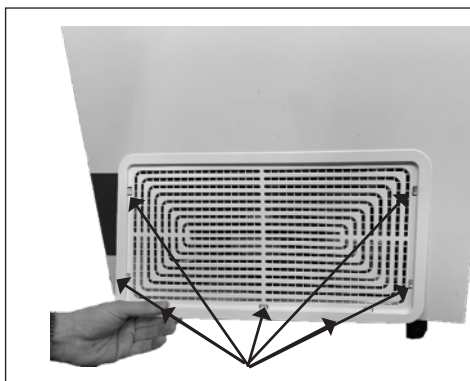2. Unlock all 7 clamps

3. Gently pull the compressor grille.


Starting device replacement

Front with terminals.

Back with connection plug.

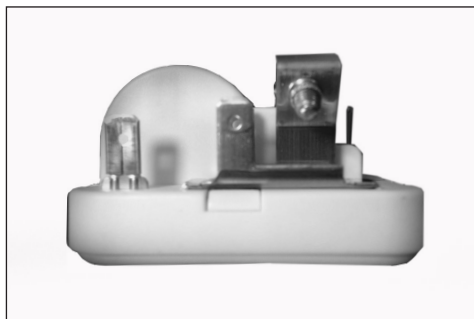Side view.

The starting device is mounted on the left side of the compressor.


1. Dismount the cover for starting device by loosen the phillips screw.

2. Push the plastic cover up.


3. Pull the cover back to loosen.

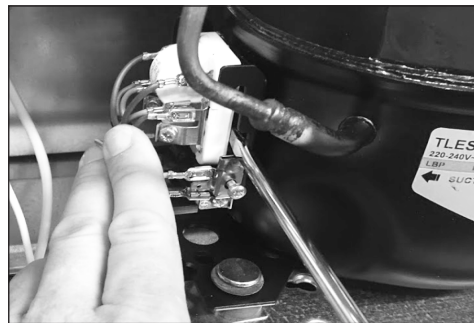
4. Use a screw driver and gently remove the starter from the socket of compressor.

5. Starting device loose from socket

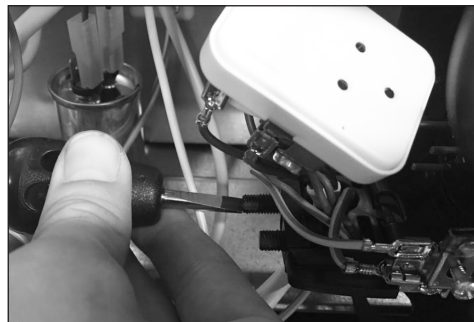6. Use a small screwdriver or a wrench size 6 – to unmount the wire fastener

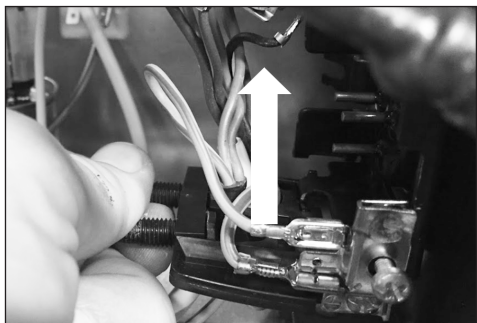7. Push the plastic bracket up.

8. Use a nose plier to unmount the wire sockets from starting device.

9. Exchange the wires 1/1 from the old starting device to the new one.

Thermostat sensor replacement


The thermostat sensor is placed inside the compartment of the appliance.


The thermostat sensor is behind the cover.


1. Take it out gently

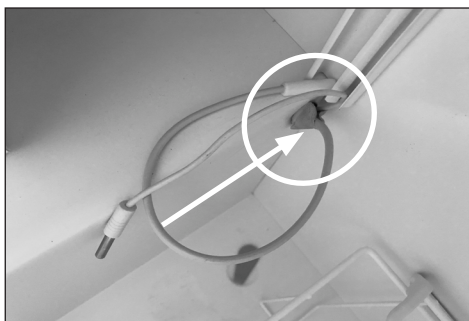
2. Remove the wire sealing and the sensor from the cover.

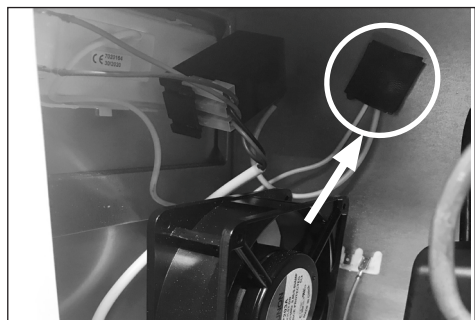
3. In Compressor compartment, remove black sealing, and gently pull the white wire until the probe is visible.

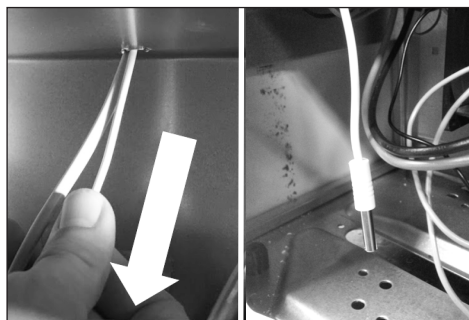
4. Gently pull the white wire until the probe is visible.

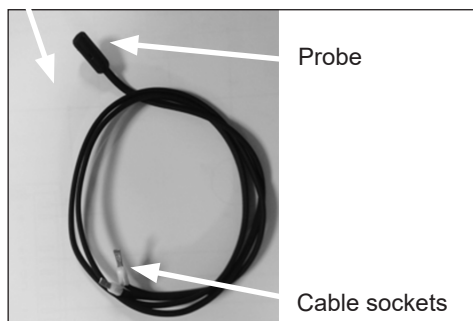5. The thermostat sensor comes with probe, wire and cable socket.

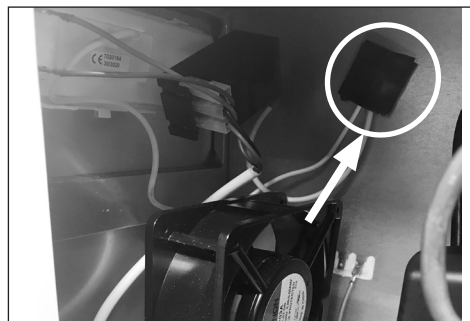6. **IMPORTANT!** When re-mounting the new thermostat sensor remember to properly seal the wire feed through.

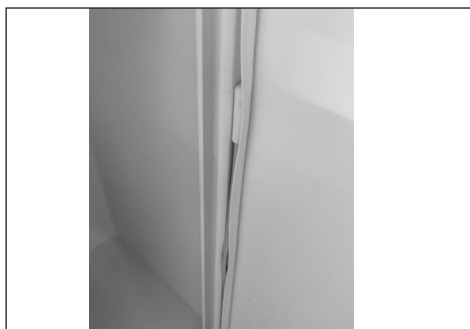7. **IMPORTANT!** When re-mounting the new thermometer make sure the wire sealing plug is placed properly.

Thermostat replacement

Control board

Use a screwdriver to gently remove the plastic cover dial.


Use a socket wrench
Size 13mm to unscrew the fixing nut


The thermostat is placed in the left corner of the compressor compartment on the terminal bracket

The thermostat is to the right of the thermometer

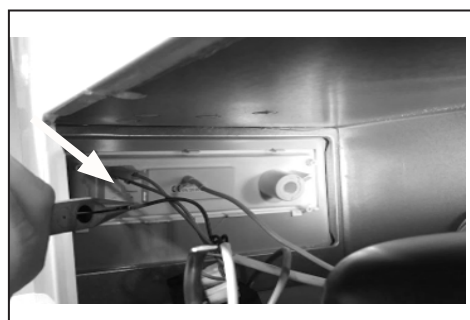Remove the thermostat from the bracket inside the compressor compartment

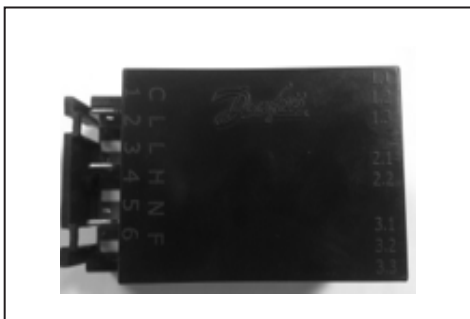Front with wire 3x sockets


Front with wire 3x sockets

Side view

Back

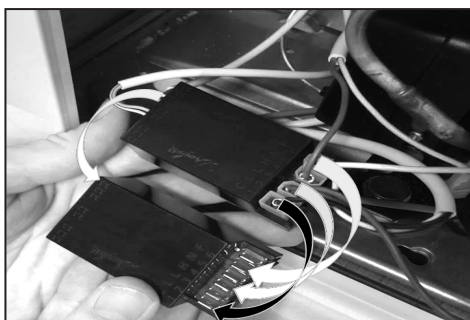Exchange the cable sockets/wires 1/1 from the old thermostat to the new

Thermometer replacement

The thermometer display is placed at the front of the appliance

Thermometer display


The thermometer sensor is placed inside appliance on the front inner-lining

Thermometer sensor is fitted into a white cable holder


Sensor probe / Grab the wire and pull it out gently

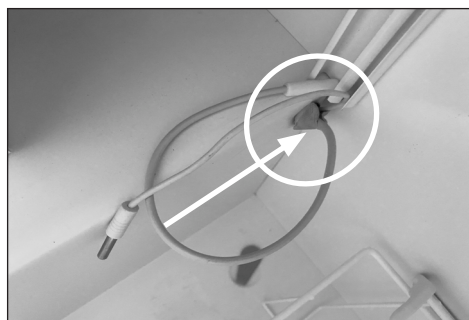2. Remove the wire sealing and the sensor from the cover.


3. In Compressor compartment, remove black sealing, and gently pull the white wire until the probe is visible.

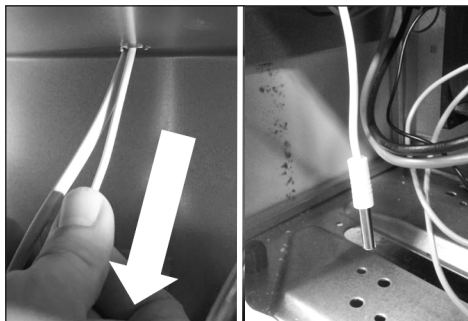4. Gently pull the white wire until the probe is visible.

Click out the thermometer display from the bracket, and pull it out with the wire

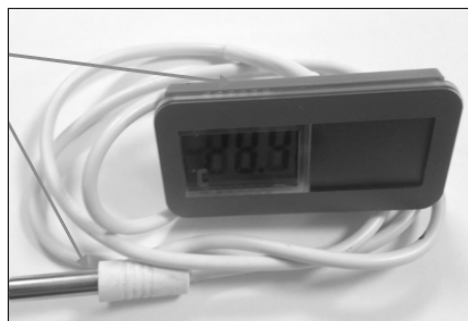
Thermometer comes with display, wire and sensor


Click the new thermometer into the bracket

Thermometer display is in place

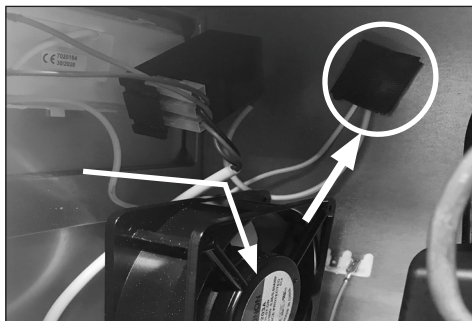IMPORTANT! When re-mounting the new thermostat sensor remember to properly seal the wire feed through.



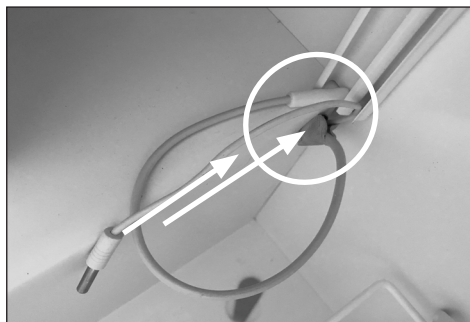
IMPORTANT! When re-mounting the new thermostat sensor remember to properly seal the wire feed through.


IMPORTANT!

When re-mounting the new thermometer make sure the wire sealing plug is placed properly.

Adjustment of lid

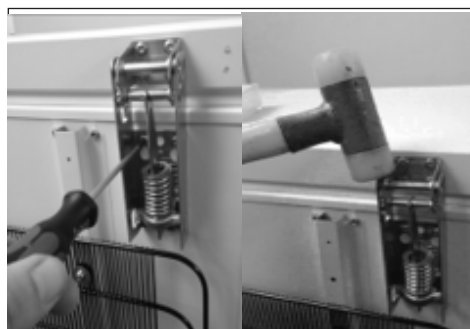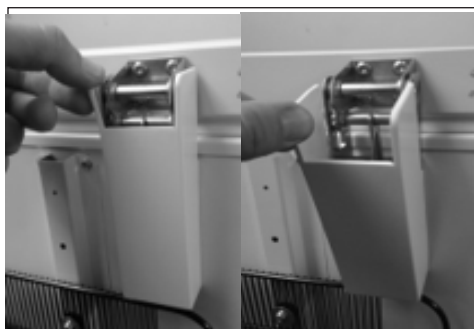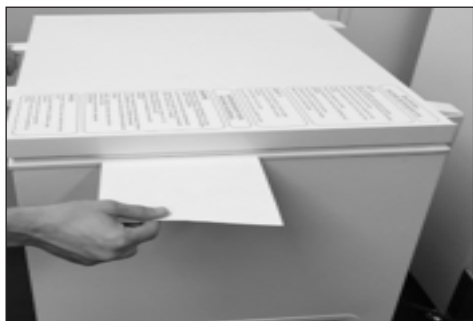Wire diagram

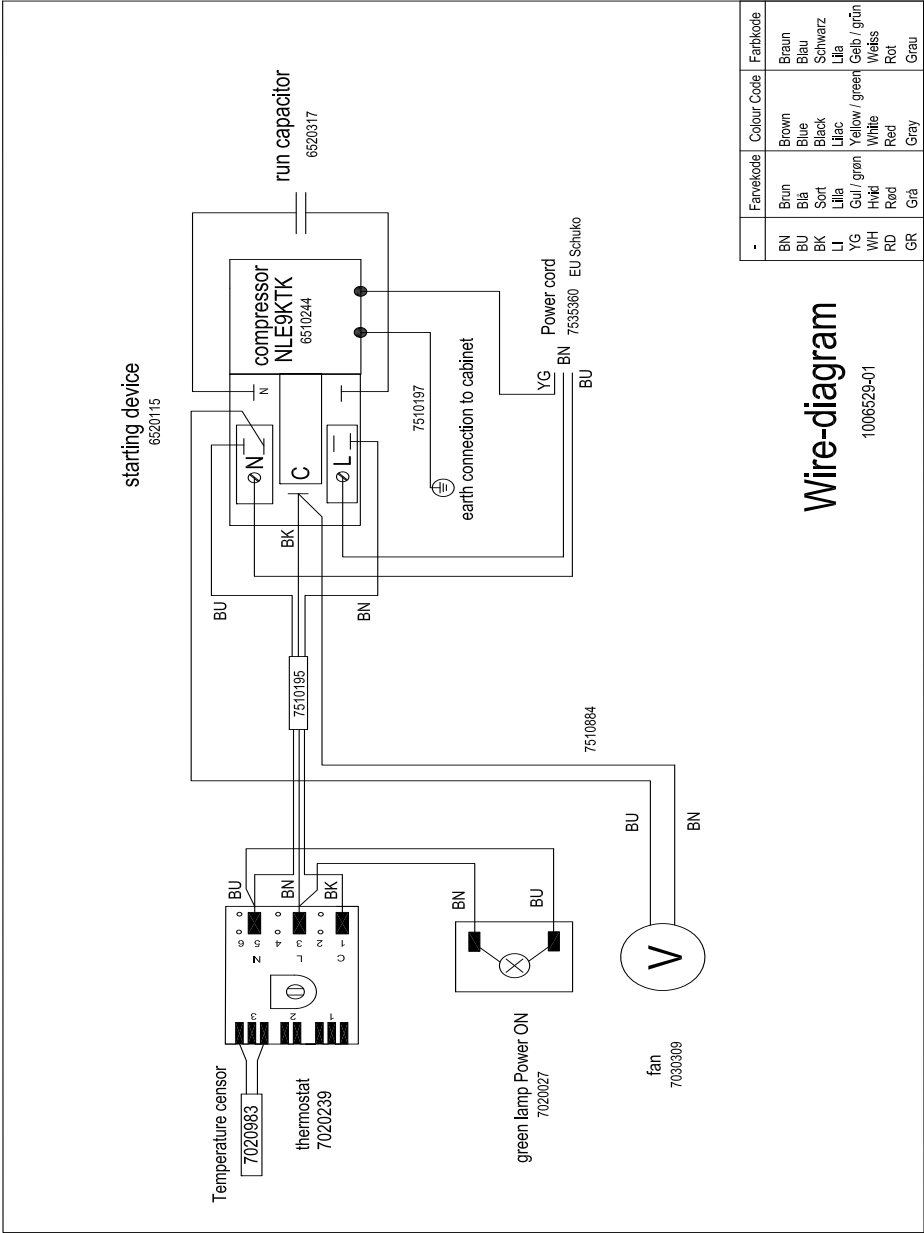Compressor replacement

Procedure of compressor switch.

1: **WARNING!** Drain coolant R600a from refrigeration system by vacuum suction.

2: **IMPORTANT!** Blow refrigeration system with NO/Nitrogen

3: Cut

A: Suction and pressure tube

B: Capillary tube

C: Dry filter

4: Dismount starting device

5: Dismount old compressor

6: Insert new compressor

7: Solder

A. Suction and pressure tube

B. Capillary tube

C. Dry filter

8: Install starting device

IMPORTANT! When solder copper tubes to iron tubes use silver tin

Filling of new refrigerant

9: Drain refrigeration system by vacuum suction

10: Check type sticker for required amount of R600a to fill on refrigerant system

On-site checklist

- Is the control panel on (Power check)
- Temperature records (manual records, FT2 data)
- Is the internal temperature inside the acceptable range of -15° to -20°
- Is the vaccine compartment clean and without ICE BUILD UP more then 0,5 cm on inner lining
- Is the Compressor running
- Is the appliance placed according to instruction in the manual.
- Does the lid close tight to cabinet and is the lid gasket in good condition
- Is the grille for compressor compartment clean
- Is all electrical components working properly
- Over all condition of the cabinet –internal and external: any corrosion, rusting, cracks
- Inspection of the refrigeration line (the condenser, evaporator, the whole refrigeration circuit/line)


Trouble shooting

Fault	Possible cause	Remedy
Compressor is not running, and the ice packs are not cold	Be patient, it is most likely that the compressor will start within a few minutes.	If this is not the case, check the following: - Check that power is connected. - Check the voltage stabilizer is properly connected and turned on. - If the above is OK, call technical supervisor.
Compressor is running, and the temperature is too high	The ventilation grille is blocked. The lid is not closed properly. The lid gasket is faulty The temperature in the room in which the appliance is installed is too high.	Ensure unhindered air circulation. Ensure that the lid is closed properly. Shield the appliance against direct sun light and ensure more ventilation to the room.
No temperature is displayed	The thermometer is broken. There is not enough light for the solar sensor.	Change the thermometer. Turn on the light.


Technical support

When contacting Vestfrost Solutions technical support please supply below information:

1. Model
2. Serial number
3. What is the issue

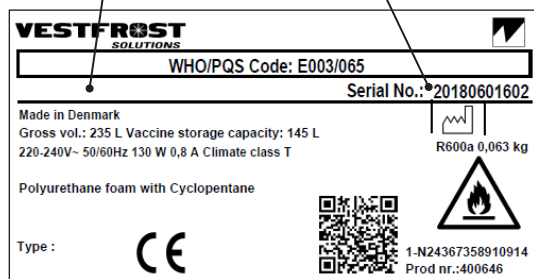Rating plate

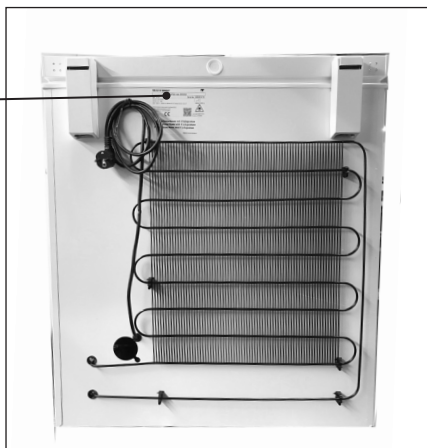Recycling procedures

Information for Users on Collection and Disposal Old Equipment and used Batteries

This symbol on the products, packaging, and/or accompanying documents mean that used electrical and electronic products and batteries should not be mixed with general household waste. For proper treatment, recovery and recycling of old products and used batteries, please take them to applicable collection points, in accordance with your national legislation and the Directives 2012/19/EU and 2006/66/EC.

By disposing of these products and batteries correctly, you will help to save valuable resources and prevent any potential negative effects on human health and the environment which could otherwise arise from inappropriate waste handling.

For more information about collection and recycling of old products and batteries, please contact your local municipality, your waste disposal service or the point of sale where you purchased the items.

Penalties may be applicable for incorrect disposal of this waste, in accordance with national legislation.

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