

ULTF-51i

ULT TECHNICAL SPECIFICATIONS

Multi voltage.



-20°C
-86°C

DIMENSIONS

Outer Dimensions HxWxD, mm	836x699x720
Depth incl. door handle	790
Inner Dimensions HxWxD, mm	310x360x445
Weight Gross/Net, kg	108 / 82
Material inner cabinet	Stainless steel
Material outer cabinet	Painted Steel
Packaging weight, kg	26
Packaging dimensions HxWxD, mm	1050x750x850
Insulation thickness	120
Insulation type	Polyurethane with Cyclopentane
Mobility	2 x adjustable feet, 2 x castors
Refrigerant, Type / gram	Nature R 2/99gr

FEATURES

Lock	☒
LED light	☐
Battery Backup for Controller, 24h	☒
Porthole	☒
Porthole size	20 mm
Dry contact	☒
Door reversibility	☐
Vacuum valve	☒
VIP (Vacuum Insulated Panel)	☒

STORAGE

Capacity (L)	50
--------------	----

ALARMS

High / Low temperature	☒
Open door	☒
Power failure	☒
Probe failure	☒

COOLING

Refrigerant	Nature R
Noise	< 48
Voltage	100-240 V
Frequency	50/60 Hz
Temperature range	-20 / -86
Ambient condition	25 °C

Products are subject to change due to enhancements and continuous development. Vestfrost Solutions reserve the right to alter any information, without further notice.

ULTF-51i

ULT TECHNICAL SPECIFICATIONS

Multi voltage.

ULTF-51i		
Voltage/Frequency	Voltage/Hz	100-240 V, 50/60Hz
Max Ambient	°C	32
Max Humidity	% rh	65
PERFORMANCE		
All data in RT20°C		
Temperature Range	°C	-20 to -86
Uniformity in performance - difference between top and bottom	°C	+/- 2,5
Pull down time	Minutes	90 min. to -75°C
Hold over time	Minutes	N/A
Noise	dB	< 48
Energy Consumption, kWh / 24h	kWh/24h	4,34 kWh/24h Set -80
Energy year	kWh/year	1584kWh/y Set -80°C
Instant Power Consumption	kW	PD 0,600/Stable 0,280
Heat Rejection*	W	200
U-Value	W/m ² K	N/A
COOLING COMPONENTS		
Refrigerant/Amount (gram)		Nature R 2/99gr
Number of compressors	pcs	1
Variable speed compressor	Yes/No	Yes
Internal air distribution (type of)		Static
Evaporator Fan	Yes/No/Variable	No
Condenser Fan	Yes/No/Variable	Yes
Number of probes	pcs	2
Defrost	Yes/No	No
CONTROLLER		
Controller	y/n/optional	Yes
Controller type		i-Care touch
USB Connection	y/n	Yes
Data connection	USB	Yes
Controller abilities		Logging of data & alarms
Controller languages		EN, DE, FR
Temperature curve in controller	y/n	Yes

* Heat Rejection is defined as average power based on energy consumption, rounded up to nearest 50W.


Products are subject to change due to enhancements and continuous development. Vestfrost Solutions reserve the right to alter any information, without further notice.

ULTF-51i

ULT TECHNICAL SPECIFICATIONS

Multi voltage.

SENSOR POSITION

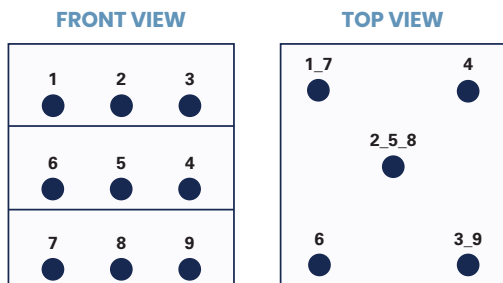

MODEL


ULTF-51i

Test type	8-point test
Test environment	Controlled conditions, empty cabinet
Ambient temperature	25°C
Humidity	60%
Set-point	-82°C
Sensor used	25gr tinned brass formed as a cylinder with a diameter of 15,2 mm
Installation	Appliance installed according to instruction manual conditions
Refrigerant	Nature R


SENSOR TEMPERATURE

Sensor position	P1	P2	P3	P4	P5	P6	P7	P8	P9								
Max	-84,2	-84,8	-85,0	-82,8	-83,3	-83,4	-80,3	-80,5	-81,0								
Avg.	-83,3	-83,8	-84,0	-82,0	-82,5	-82,5	-79,4	-79,5	-80,0								
Min.	-83,7	-84,2	-84,4	-82,3	-82,8	-82,8	-79,7	-79,9	-80,4								

WARM UP & PULL DOWN

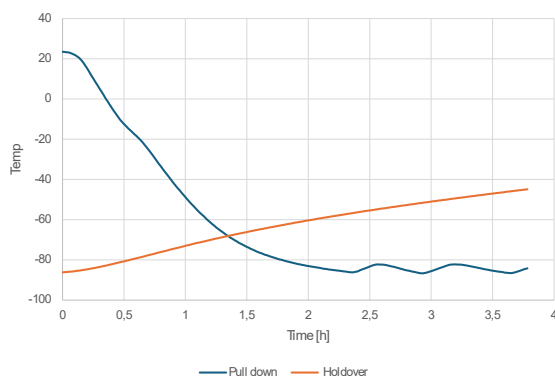

CYCLIC OPERATION

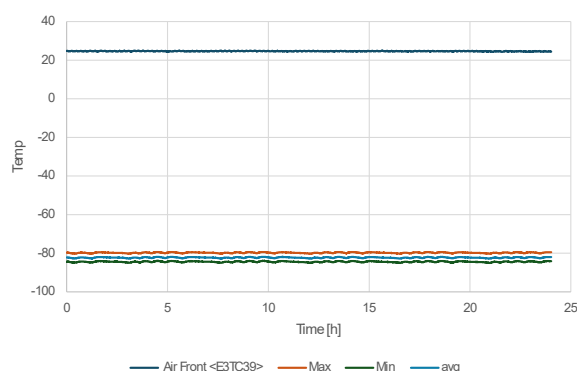

TYPICAL PERFORMANCE IN AMBIENT 25°C - EMPTY CABINET

Avg. cabinet temperature at -82°C setpoint	-82,8°C
Uniformity	+/- 2,5°C
Stability	1,1°C
1 min. door open recovery to -75°C	32 min.
Cycle rate -82°C setpoint on/off	50/5 min./min.
Duty cycle -82°C setpoint	90%
Energy consumption	4,35 kWh/day
Pull down time to -75°C avg. temperature	103 min.
Hold over time from avg. -82°C to -60°C	102 min.
Heat rejection	200 W

Products are subject to change due to enhancements and continuous development. Vestfrost Solutions reserve the right to alter any information, without further notice.